


Medical Mask



Description	<p>This medical mask is made up of three skin-friendly layers which protect against droplets that could potentially carry harmful particles and bacteria. Thanks to the flexible material at the top of the mask, it is easily adjustable to the wearer's face and can be comfortably worn with glasses and goggles. Despite multi-layer processing, the mask is characterized by high air permeability.</p> <p>Manufactured in the EU.</p> <p>Certified for use in the medical field.</p>
Classification	Type II
Usage	Type II is suitable for use in the medical and surgical field. It offers high antibacterial and antiviral protection provided that basic hygiene rules are maintained and the mask is replaced regularly after a certain period of time.
Material	<ol style="list-style-type: none"> 1. Outer layer: spunbonded propylene nonwoven Performs a protective role and stops small impurities. 2. Middle layer: meltblown filtering nonwoven Blocks dust particles and aerosol-borne bacteria from entering, has a high resistance to water pressure, good air permeability and a good resistance to acids and alkalis. Suitable for high and low temperatures. 3. Inner layer: spunbond propylene nonwoven Provides good ventilation that allows free breathing. 4. Nose wire Provides flexible adjustment to the shape of the face.
Bacterial filtration efficiency (BFE)	$\geq 99\%$
Breathability	$\Delta P < 40 \text{ Pa/cm}^2$
Certification	DIN EN 14683 + AC: 2019-09 CE-certified

Product Characteristics

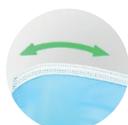


The nose wire at the top of the mask is made of flexible material. It increases comfort and prevents glasses from fogging.



Elastic and soft straps ensure high comfort, even if worn for long periods of time. Straps sewn into the mask keep it firmly on the face, preventing slipping.



Despite multi-layer processing, the mask is characterized by high air permeability, enabling free breathing.

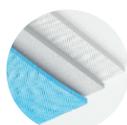

The mask is made up of different polypropylene layers with filters that effectively trap particles in the air. The mask is odorless, highly air-permeable and does not impede breathing.

How to use

1.



Open the mask and pull the ear straps apart with both hands.

2.



Hold the mask against your chin to completely cover your nose and mouth.

3.



Pull the ear straps behind your ears and adjust them for comfortability.

4.



Use both hands to press on the middle of the nose wire. Press the flexible material onto the face, moving inward and on both sides of the nose, until it sits comfortably and firmly.

To place an order, please get in touch via email: sampling@inveox.com

Due to high demand and the need to integrate and optimize various supply chains, the delivered product may vary in appearance and packaging from the item in this description.

inveox GmbH • Lichtenbergstraße 8, 85748 Garching bei München, Germany
www.inveox.com • Phone: +49 (0) 89 / 57 84 76 01 • Fax: +49 (0) 89 / 41 41 60 543